


Bird flu, Avian influenza Rapid Test Kit



15min Rapid Test Kit to detect avian influenza antigen in swabs of throat, feces, etc.

Ringbio® Avian influenza Antigen Rapid Test Kit is based on lateral flow immunoassay to detect nucleoprotein (NP) of AIV as target in avian samples. The total detection time is 5-15min. It is suitable for rapid detection of bird flu in various samples of wild and domestic birds.

Kit Components

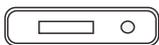

Test Card, 25 pcs



Pipette & swab, 25 pcs



Buffer tube, 25 pcs



Instruction manual, 1 set

Operations

15min operation

1




Take Sample

2



Extract with buffer

3

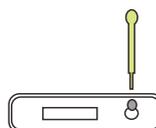


Add to test card

4




Wait for 15min

5

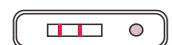

Check Result

**Collect samples with the swab provided in the kit.

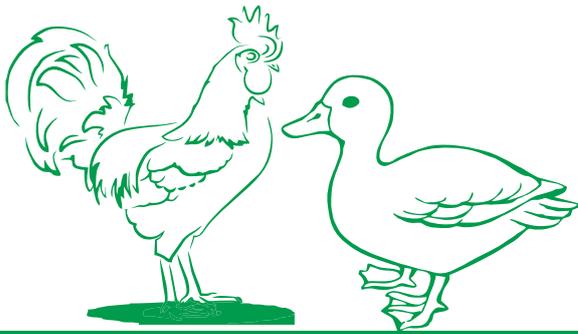

Avian Influenza Ag Rapid Test Kit

Technical Specifications

Product Number	T21195
Unit Size	25 tests / kit, 50 tests / kit
Principle	Lateral flow immunoassay
Applied Sample	Throat swabs
Species	Avian, poultry, eagle, owl, and other wild birds
Storage	2-30 °C
Shelf life	12 months
Specificity	>95% compared with ELISA, RT-PCR
Sensitivity	>90% compared with HI/HA
Shipment	Room temperature

Other Influenza Test Kits



AIV H5 Antibody ELISA Kit

Product Number: E20193, 96 tests, 480 tests
Application: poultry serum



Bird flu, Avian influenza Real-time qPCR Kit

Product Number: QP008, 50 tests
Application: swabs, tissue, feces.

You can always visit www.ringbio.com for our latest product update
Not all test kits are registered in all countries, consult your sales before order
For veterinary use only, for research use only, for outside of China only.